

# LetsGetChecked Blood Home Sample Collection IFU 203 SPECS AND INSTRUCTIONS

## Colour Reference:

■ LGC Teal: Pantone 7710 CP (CMYK)

■ LGC Navy: Pantone 654 CP (CMYK)

■ LGC Red: Pantone 1925 CP (CMYK)

■ CMYK

## Size:

108.5mm wide by 152.5mm high **folded**

## Paper:

135gsm silk

## Fold type:

5 page accordion

## Font:

Proxima Nova

(Regular, Semibold and Bold)

## Part Number:

QR-IFU-203-R1 Jun-2025




\*Visual reference for how to fold down the sheet so the front page is the first one to see. Artwork and page count may vary.


## FOLD GUIDE

----- Fold


## Accordion fold

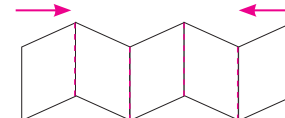

## Inside

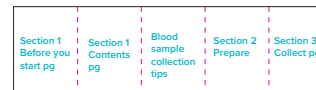

## Outside

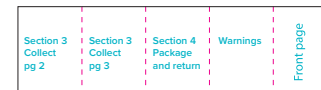

# 1 Before you start



Some laboratory tests require samples collected while fasting or at a specific time. Visit [LetsGetChecked.com/sample-collection](https://letsgetchecked.com/sample-collection) to learn what you need to do for your test.



Collect your sample in the morning Monday-Friday and return on the same day.



Plan your sample return in advance. Check the return delivery card in your kit for details.



Download the LetsGetChecked App from your app store.



Scan the QR code below for an instructional video.



## Check your kit contents\*



Gauze pack



Alcohol swabs × 2



Bandages × 2



Blood collection tube




Sterile lancets × 4



Biohazard bag



Return envelope

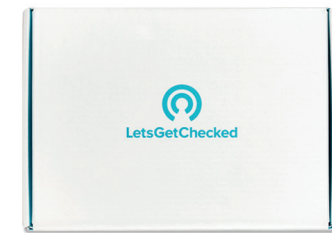

Kit box

\*The look or color of items may vary from what you see in these images.

## Blood sample collection tips



**Collect from warm hands**  
Run your hands under warm water.



**Get the blood flowing**  
Light exercise such as walking can help.



**Stay hydrated**  
Drink water before you start.



**Check your nails**  
Collection may be easier with short nails.



**Keep your arm angled downwards**  
Keep your elbow below your heart and your hand at waist level or lower.



**Select your finger**  
We recommend using your 4th (ring) finger on your non-dominant hand. Avoid fingertips with hard skin or calluses. Do not pierce the baby/pinky finger.

# 2 Prepare for sample collection



Activate your kit by scanning the QR code on the biohazard bag or on the **LetsGetChecked App** or [LetsGetChecked.com/activate](https://letsgetchecked.com/activate). Only activate when you are ready to collect your sample.



Wash and dry your hands.



Remove all the contents from your kit box. Open the packaging.

# 3 Collect your blood sample



Read the blood sample collection tips. Twist the cap to open the blood collection tube.



Twist off the white lancet cap. Do not pull it.



Clean your finger with the alcohol swab and allow to air-dry.

Please continue overleaf.





When you're ready, position the lancet to one side of the fingertip. Press down firmly until you hear a click.



Use the gauze to wipe away the first small drop of blood.



Massage your hand and finger. Do not touch or squeeze the fingertip — this prevents blood flow.



Fill the blood collection tube to the top line. An under-filled tube will not be accepted.




Apply a bandage to your finger.



Replace the cap on the blood collection tube. Press down firmly until you hear a click.




Lightly tap the blood collection tube with your finger. Turn the tube upside down ten times to mix the sample. This preserves your sample.

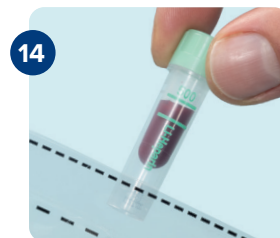

Place the blood collection tube in the biohazard bag and seal the bag.

## 4 Package and return

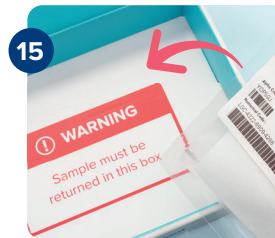

Place the biohazard bag into the kit box. Double check you have activated your kit.



Place the kit box in the return envelope provided.



Send the kit back on the same day. Check the return delivery card in your kit box for your return method.

### ⚠ WARNING:

- You must be over the age of 18 to use this home sample collection kit.
  - The sample and information provided must be your own. Test results may not be released to a third party.
  - Failure to follow these instructions may impact the accuracy of your test result or cause your sample to be rejected by the laboratory.
  - Some medication and supplements may impact your test result. Do not make any changes to your medication without consulting your healthcare provider.
  - If you are collecting a sample for a test that requires fasting please consult with your healthcare provider, particularly if you have any medical conditions.
  - Do not pierce or collect blood from the smallest finger or any other body part not indicated in these instructions.
  - Do not use lancets on damaged, scarred, or infected skin.
  - Lancets may cause some bruising and tenderness at the puncture site. There is a small possibility of pain, nerve damage, scarring, bone damage, or infection. Consult your healthcare provider if you are concerned.
  - There is a risk of fainting with any blood sample. If you are concerned, consult your healthcare provider before collecting your sample. If you proceed, ensure you are seated and have someone with you. If you feel faint or unwell during sample collection, stop and seek medical advice.
  - Use only the contents included in this single-use kit for collecting your sample. Check the expiry date before you use the kit and do not use if expired.
  - If the kit or any of the contents are damaged or missing, the kit should not be used. Do not tamper with the kit contents — risk of serious injury.
  - Used kit contents can generally be safely disposed of in household waste. Refer to your local laws and guidelines.
  - Choking hazard — this kit contains small parts which may present a choking hazard. Keep away from children.
  - The sample must be returned in the kit box.
  - If you have difficulty using this kit due to vision impairment or difficulties in handling the test components, or understanding test instructions or results, please contact a healthcare provider.
  - This test system is not a substitute for visits to a healthcare provider. The information provided by this product should not be used to start, stop, or change any course of treatment unless advised by your healthcare provider.
- Do not use this kit if:**
- You have a history of fainting or a medical condition that increases your risk of fainting.
  - You have a bleeding disorder or are taking medication that increases the risk of bleeding.
  - You have impaired lymphatic drainage in both upper limbs.
  - You have a condition associated with poor healing, increased risk of infection, or skin ulceration — including but not limited to Raynaud's phenomenon, scleroderma, peripheral vascular disease, or peripheral neuropathy.
- Purpose/Intended use of device:**  
The LetsGetChecked Blood Home Sample Collection Kit is an in vitro diagnostic medical device specifically intended for the primary containment and preservation of blood samples, for the purpose of in vitro diagnostic examination by our partner accredited laboratories.
- Contact Customer Care:**  
1-855-554-1959 (TTY 711)

PrivaPath Diagnostics Ltd.,  
Unit 1, Northern Cross Business Park,  
North Road, Dublin 11,  
D11 XT26, Ireland.

Storage temperature: 15°C — 30°C.  
Transport temperature: 4°C — 25°C.



For explanation of symbols see  
[LetsGetChecked.com/symbol-glossary](https://www.letsgetchecked.com/symbol-glossary)

LetsGet  
Checked

LETSGETCHECKED  
Blood Home Sample  
Collection Kit

# How to collect your sample



⚠ Activate your kit before you collect your sample