LetsGetChecked Optum House Calls
Stool Home Sample Collection IFU 242
SPECS AND INSTRUCTIONS

Colour Reference:
■ LGC Teal: Pantone 7710 CP (CMYK)
■ LGC Navy: Pantone 654 CP (CMYK)
■ LGC Red: Pantone 1925 CP (CMYK)
■ CMYK

Size:
108.5mm wide by 152.5mm high **folded**

Paper:
135gsm silk

Fold type:
4 page accordion

Font:
Proxima Nova
(Regular, Semibold and Bold)

Part Number:
QR-IFU-242-R1 Feb-2025


*Visual reference (showing the English side of the IFU) for how to fold down the sheet so the front page is the first one to see. Artwork and page count may vary.

FOLD GUIDE



----- Fold

Fold in a half first

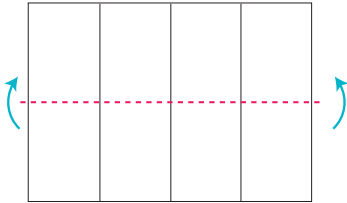Then accordion fold with the text in English on the outside

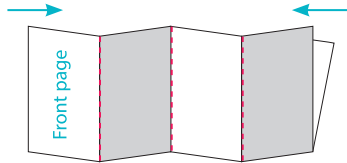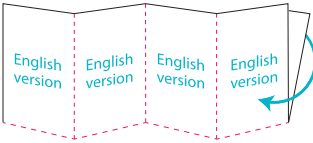Outside – English version

EN	EN	EN	EN
Front page	Section 1 Before you start pg	Section 1 Contents pg	Section 2 Prepare
EN	EN	EN	EN
Section 3 Collect	Section 4 Package and return pg 1	Section 4 Package and return pg 2	Warnings

Inside – Spanish version

ES	ES	ES	ES
Front page	Section 1 Before you start pg	Section 1 Contents pg	Section 2 Prepare
ES	ES	ES	ES
Section 3 Collect	Section 4 Package and return pg 1	Section 4 Package and return pg 2	Warnings

LETSGETCHECKED
Stool Home Sample
Collection Kit

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

QR-IFU-242-R1 Feb-2025

1 Before you start

-  Do not collect a sample if you are on your period or bleeding due to hemorrhoids or constipation.
-  Collect your sample in the morning and return on the same day.
-  Plan your sample return in advance. Check the return delivery card in your kit for details.
-  Scan the QR code below for an instructional video.

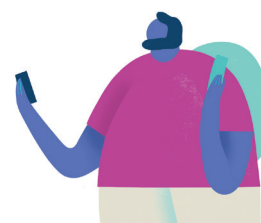Check your kit contents*

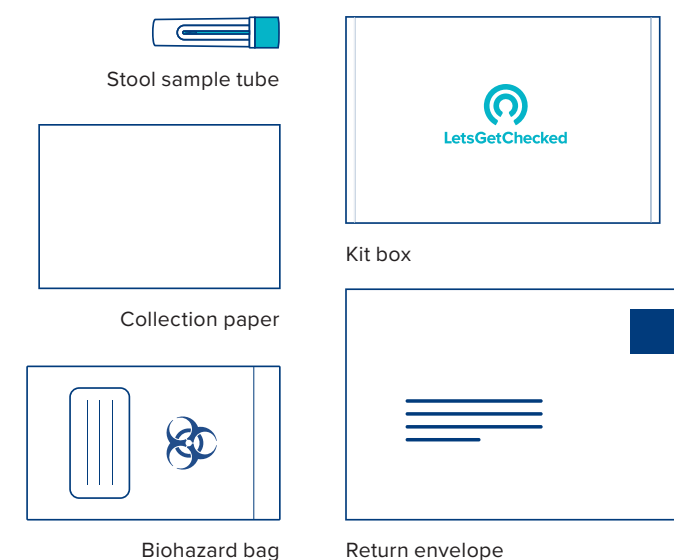
*The look or color of items may vary from what you see in these images.

2 Prepare for sample collection

Fill in the label on the biohazard bag.

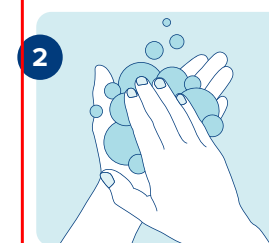
Wash and dry your hands.

TIP:


Place the kit in your restroom so it is at hand when you need it.

↓ Flip over for next steps. ↓


3 Collect your stool sample

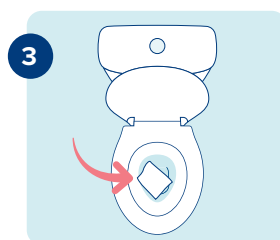Place the collection paper inside the toilet bowl and on top of the water. Pass your stool (poop) onto the collection paper.

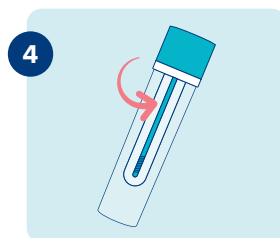
Collect your sample before the paper sinks and the stool touches water. Twist the cap to open the stool sample tube and remove the stick.

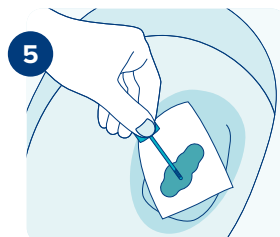Scrape the surface of the stool with the stick. Cover the grooved portion of the stick with stool sample.

4 Package and return

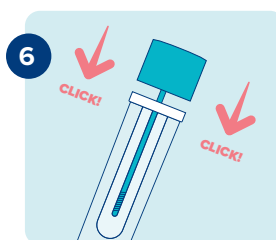Replace the stick in the stool sample tube. Press down firmly on the cap until you hear a click to seal it shut. Do not reopen the tube.

Place the stool sample tube in the biohazard bag and seal the bag.

Place the biohazard bag into the kit box. Double check you have filled in the biohazard bag label.

Flush the toilet. The collection paper is flushable.

Wash and dry your hands.

Place the kit box in the return envelope provided.

Send the kit back on the same day. Check the return delivery card in your kit box for your return method.

- WARNING:**
- You must be over the age of 18 to use this home sample collection kit.
 - The sample and information provided must be your own. Test results may not be released to a third party.
 - Failure to follow these instructions may impact the accuracy of your test result or cause your sample to be rejected by the laboratory.
 - Urine and excessive dilution of samples with water from the toilet bowl may cause erroneous test results. For best results, use the collection paper provided.
 - This kit allows for the collection of a stool sample for the assessment of human blood in the stool in our partner laboratories. Blood in the stool can happen for many reasons, one of which is colorectal cancer.
 - The test performed by our partner laboratory cannot tell if you have cancer and is not a replacement for a colonoscopy. There is a possibility that your result may be negative even if you have colon cancer. Please talk to your healthcare provider if you have any symptoms or if you are concerned.
 - Do not collect a sample if menstruating or if actively bleeding for example, as a result of hemorrhoids, constipation or urinary bleeding.
 - Do not remove liquid from the stool sample tube at any point. Do not drink or breathe in the liquid or let the liquid touch your skin or eyes. If you breathe in the liquid get fresh air straight away. If any exposure to the liquid

occurs, wash the area immediately with lots of water and consult your healthcare provider.

- Use only the contents included in this single-use kit for collecting your sample. Check the expiry date before you use the kit and do not use if expired.
- If the kit or any of the components are damaged or missing, the kit should not be used.
- Choking hazard – this kit contains small parts which may present a choking hazard. Keep away from children.
- The sample must be returned in the kit box.
- If you have difficulty using this kit due to vision impairment or difficulties in handling the test components, or understanding test instructions or results, please contact a healthcare provider.
- This test system is not a substitute for visits to a healthcare provider. The information provided by this product should not be used to start, stop, or change any course of treatment unless advised by your healthcare provider.

Purpose/Intended use of device:

The LetsGetChecked Stool Home Sample Collection Kit is an in vitro diagnostic medical device specifically intended for the primary containment and preservation of stool samples, for the purpose of in vitro diagnostic examination by our partner accredited laboratories.

Contact Customer Care:
1-855-554-1959 (TTY 711)

 PrivaPath Diagnostics Ltd.,
Unit 1, Northern Cross Business Park,
North Road, Dublin 11,
D11 XT26, Ireland.

 Storage temperature: 15°C — 25°C.
 Transport temperature: 1°C — 30°C.

For explanation of symbols see
[LetsGetChecked.com/symbol-glossary](https://www.letsgetchecked.com/symbol-glossary)

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