

LetsGetChecked Blood Sample ImPress IFU for
demo use - QR-IFU-123-R0

SPECS AND INSTRUCTIONS

Colour Reference:
CMYK

Size:
108.5mm wide by 152.5mm high **folded**

Paper:
135gsm silk

Fold type:
4 page accordion

Font:
Proxima Nova
(Regular, Semibold and Bold)

Part Number:
QR-IFU-123-R0 Jun-2025




*Visual reference for how to fold down the sheet so the front page is the first one to see. Artwork and page count may vary.

FOLD GUIDE

----- Fold

Accordion fold

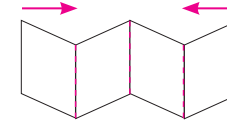

Inside







Section 1 Before you start pg	Section 1 Contents pg	Section 2 Collect sample steps pg 1	Section 2 Collect sample steps pg 2
-------------------------------------	-----------------------------	--	--


Outside

Section 2 Collect sample steps pg 3	Section 2 Collect sample steps pg 4	Warnings	Front page
--	--	----------	------------


1 Before you start

Estimated time needed: 10-15 minutes.

-  Read all the instructions including the warnings.
-  An instructional video can be viewed at [LetsGetChecked.com/instructions-imp-collection](https://letsgetchecked.com/instructions-imp-collection)
-  This is a single use device. Once you have removed the tabs, do not press the device until it is positioned on the arm and you are ready to begin sample collection.
-  Short sleeves are recommended for sample collection.
-  Have a timer ready to set for 10 minutes.
-  Contact customer support at [help.LetsGetChecked.com](https://help.letsgetchecked.com) if the device does not function as described.

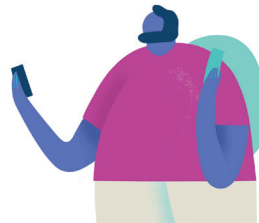

Check your kit contents

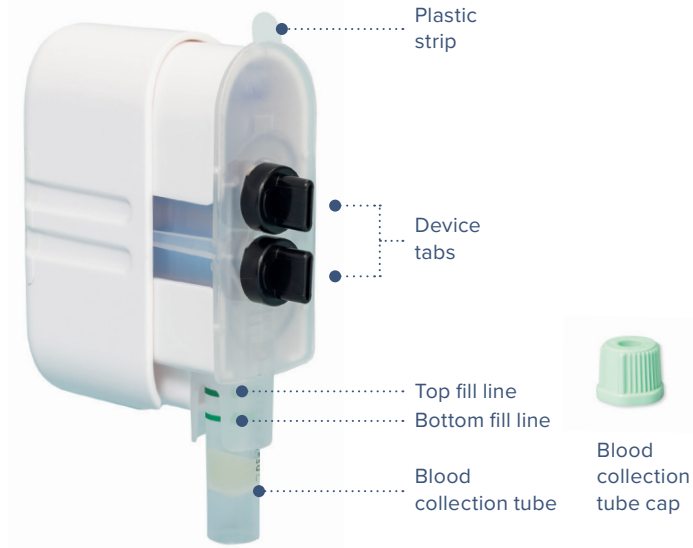

Also required but not provided:

- ☐ Timer
- ☐ Bandages
- ☐ Sterile Gauze
- ☐ Alcohol swab
- ☐ Medical gloves (if collecting the sample from someone else)

2 Collect your blood sample




Wash and dry your hands. If you are collecting the sample from someone else, the use of medical gloves is recommended.




Remove the device and tube cap from the packaging. You will use the cap later to seal the sample in the blood collection tube.



Clean the upper arm with an alcohol swab and allow to air-dry.

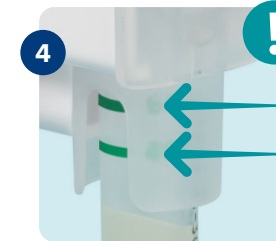


Take note of the two fill lines on the tube. You will aim to fill to the top line with blood. The top fill line may be hard to see once the device is placed on your arm.

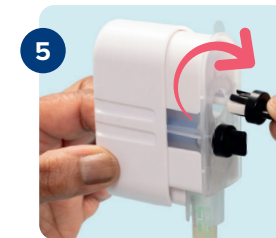

Twist and remove the black tabs from the device.



Peel off the clear plastic strip using the tab at the top of the device.

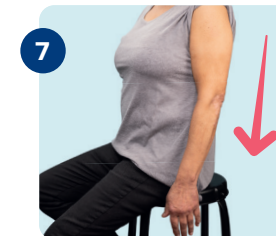

Keep the arm straight and fingers pointing down towards the floor. Perform blood collection in an upright seated position.

Please continue overleaf.




Stick the device on the upper arm with the collection tube pointing downwards.



Using your palm, firmly press the device against the arm until you hear a click. Release to start the collection process.



Start a timer for 10 minutes. It can take up to two minutes to see the first drop of blood in the collection tube.

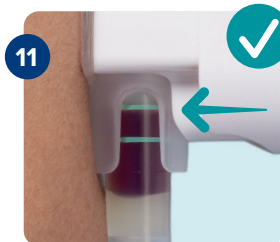

Follow the next step to remove the device when:

- Blood has filled to the top line or
- 10 minutes have passed.



Peel the device from the arm while keeping the tube upright.



To remove the tube from the device, gently pull it downwards. Discard the device.



Place the cap on the blood collection tube and press down firmly until you hear a click. Clean arm with gauze and apply a bandage if needed.



Turn the tube upside down 10 times to mix the sample. Gently tap or flick the tube to get blood moving if needed. Wash and dry your hands.

⚠ WARNING:

- You must be over the age of 18 to use this device to self-collect a blood sample.
- Single use only. Not intended for more than one use. For use only on a single patient.
- Not for use on infant heels.
- Do not use after the expiration date.
- Do not resterilize.
- For external use only.
- Use personal protective equipment when handling biohazardous materials.
- The lancet is single use and automatically retracts into the device in a safe position after use.
- Improper use of blood lancets can increase the risk of inadvertent transmission of bloodborne pathogens, particularly in settings where multiple patients are tested.
- Minor bruising or residual marks may occur at the sample collection site. Consult with a healthcare provider if you are concerned or if this persists.
- Do not use the device on damaged, scarred, or infected skin. There is a small possibility of pain, nerve damage, scarring, bone damage, or infection. Consult a healthcare provider if you are concerned.
- Use while seated as fainting may occur with any blood sampling procedure. If you are concerned, consult a healthcare provider before collecting your sample. If you proceed, ensure you have someone with you. If you feel faint or unwell during sample collection, stop sample collection. Consult a healthcare provider if you continue to feel faint or unwell.
- Do not use if any of the contents are damaged or missing. Do not use if device packaging has been opened or damaged. Do not tamper with the device — risk of serious injury.
- Discard the entire device after use. Follow local guidelines for discarding used sharps as a biohazard material.
- Laboratories must validate use of product for their specific diagnostic assay.
- Choking hazard — keep out of reach of children.
- Contact our team at +353 1 513 6654 before using this device if you have any questions or concerns.

Do not use this device without consulting a healthcare provider if:

- You have a history of fainting or a medical condition that increases your risk of fainting.
- You have a bleeding disorder.
- You have a condition associated with an increased risk of infection, impaired wound healing, skin ulceration, or swelling in the arms.

Purpose/Intended use of device:

The ImPress device is intended for the collection of capillary blood samples from the upper arm into a compatible sample collection reservoir, for use by healthcare professionals or lay users aged 18 years and older.

PrivaPath Diagnostics Ltd.,
Unit 1, Northern Cross Business Park,
North Road, Dublin 11,
D11 XT26, Ireland.

Temperature limit: 15°C — 25°C.



For explanation of symbols see
[LetsGetChecked.com/symbol-glossary](https://www.letsgetchecked.com/symbol-glossary)

**LetsGet
Checked**

ImPress Blood Collection

For demonstration purposes only -
not for human use



QR-IFU-123-R0 Jun-2025