

LetsGetChecked myGene DNA Saliva Collection  
Activation and Return Guide  
SPECS AND INSTRUCTIONS

**Colour Reference:**  
CMYK

**Size:**  
108.5mm wide by 152.5mm high **folded**

**Paper:**  
135gsm silk

**Fold type:**  
3 page accordion

**Font:**  
Proxima Nova  
(Regular, Semibold and Bold)

**Part Number:**  
QR-IFU-115-R0 May-2024




\*Visual reference for how to fold down the sheet so the front page is the first one to see. Artwork and page count may vary.



FOLD GUIDE

----- Fold

Accordion fold

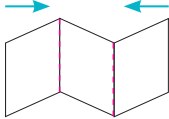

Inside



Outside




# ! Before you start

- ☐ Do not eat, drink, smoke, or chew gum for 30 minutes before collecting your sample.
- ☐ Collect your sample on a weekday and return on the same day.
- ☐ Do not remove the plastic film from the funnel lid.
- ☐ Plan your sample return in advance. Check the return delivery card in your kit for details.
- ☐ Read all instructions and warnings on the ORAgene®•Dx saliva collection kit.



## Kit contents

Make sure you have all these items in your kit:\*

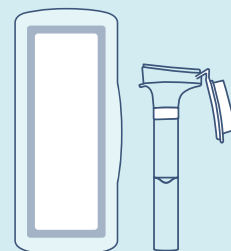


ORAgene®•Dx saliva collection kit



Biohazard bag



Return delivery card

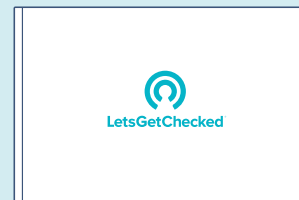

Kit box



Return envelope

\*The look or color of items may vary from what you see in these images.

# 1 Activate your kit



Activate your kit at **LetsGetChecked.com/activate**. Enter the unique alpha and numerical codes on your biohazard bag. Only activate when you are ready to collect your sample.



Fill in the label on the biohazard bag.



Wash and dry your hands.


Please continue overleaf.




## 2 Prepare for sample collection



Remove all the contents from your kit box and open the packaging.

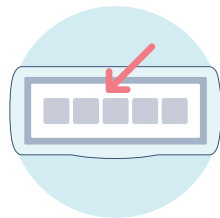

Follow the instructions inside the ORAgene®•Dx packaging to collect your sample.

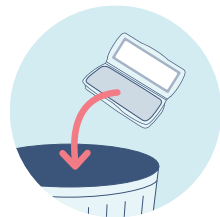

Discard the ORAgene®•Dx packaging.

## 3 Package and return



Place the sample into the biohazard bag. Then, place the biohazard bag into the kit box. Double check you have filled in the biohazard bag label.



Place the kit box into the return envelope provided.



Send back on the same day. Check the return delivery card in your kit box for your return method.



# How to activate and return your sample



⚠️ Activate your kit at [LetsGetChecked.com/activate](https://LetsGetChecked.com/activate) before collecting your sample