P1

LetsGetChecked Verily Blood DBS Home Sample Collection IFU 021

SPECS AND INSTRUCTIONS

Colour Reference:
CMYK

Size:
108.5mm wide by 152.5mm high **folded**

Paper:
115gsm silk

Fold type:
5 page accordion

Font:
Proxima Nova
(Regular, Semibold and Bold)

Part Number:
QR-IFU-021-R0 Aug-2024



*Visual reference (showing the English side of the IFU) for how to fold down the sheet so the front page is the first one to see. Artwork and page count may vary.


FOLD GUIDE

----- Fold


Accordion fold first

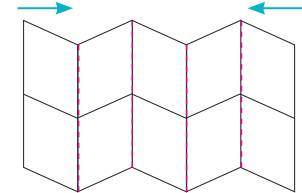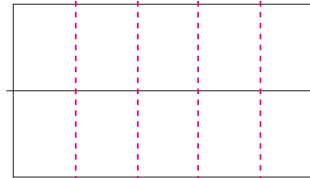
Then fold in a half with the text in English on the outside

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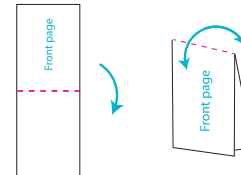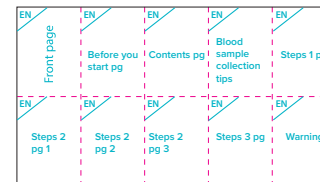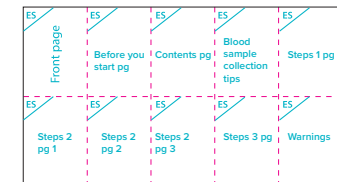QR-IFU-021-R0 Aug 2024


Before you start

- Collect the sample in the morning Monday-Friday and return on the same day.
- Plan the sample return in advance. Check the return delivery card in your kit for details.
- Read all the instructions including the warnings.

Kit contents

Make sure you have all these items in your kit.*

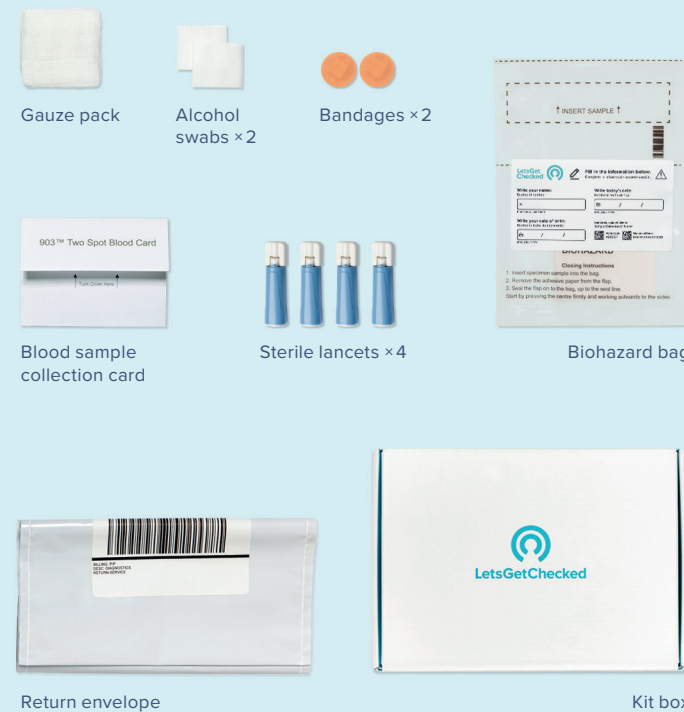*The look or color of items may vary from what you see in these images.

Blood sample collection tips

- Collect from warm hands**
Run the hands under warm water.
- Get the blood flowing**
Light exercise such as walking can help.
- Stay hydrated**
Ensure the person giving the sample drinks water before starting.
- Check the nails of the person giving the sample**
Collection may be easier with short nails.
- Keep the arm angled downwards**
The person giving the sample should keep their elbow below their heart and their hand at waist level or lower.
- Select the finger**
We recommend using the 4th (ring) finger on the non-dominant hand. Avoid fingertips with hard skin or calluses. Do not pierce the baby/pinky finger.

1 Prepare for sample collection

Fill in the label on the biohazard bag with the details of the person giving the sample.


Wash and dry the hands of the person giving the sample. If you are assisting with a child's sample collection, ensure your hands are washed also.

Remove all the contents from your kit box and open the packaging. The blood sample collection card must be kept on a clean dry surface. Do not touch the two circles.

2 Collect your blood sample

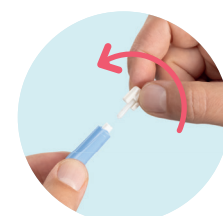
Read the blood sample collection tips. Twist off the white lancet cap. Do not pull it. Clean the finger of the person giving the sample with the alcohol swab and allow to air-dry.


When you're ready, position the lancet to one side of the fingertip. Press down firmly until you hear a click.



Use the gauze to wipe away the first small drop of blood.

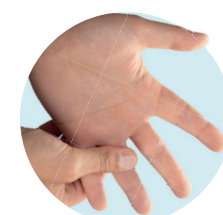Massage the hand and finger. Do not touch or squeeze the fingertip — this prevents blood flow.

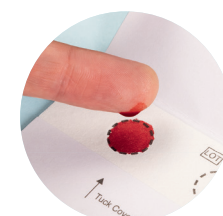Drop blood onto a circle and fill it entirely as shown. The front and back of the circle should be completely red when complete. Repeat for the second circle. At least one circle must be completely filled.

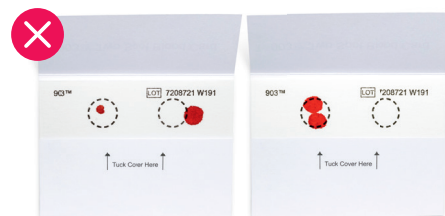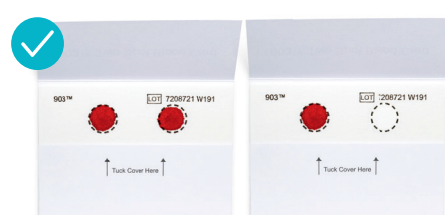
Apply a bandage to the finger.

Fold down the cover of the blood sample collection card and gently tuck it under the flap (marked with arrows). Do not flatten or press down on the card.

Allow the blood sample collection card to dry for 3 hours. Do not apply heat to assist drying.

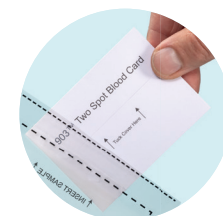Place the blood sample collection card in the biohazard bag and seal the bag.

3 Package and return

Place the biohazard bag into the kit box. Double check you have filled in the biohazard bag label.

Place the kit box in the return envelope provided.

Send the kit back on the same day. Check the return delivery card in your kit box for your return method.

WARNING:

- You must be over the age of 18 to use this home sample collection kit for yourself or to assist the collection of a child's sample.
 - The sample and information provided must be your own if you are taking the test or your child's own if they are taking the test. Test results may not be released to a third party.
 - Failure to follow these instructions may impact the accuracy of the test result or cause the sample to be rejected by the laboratory.
 - Some medication and supplements may impact the test result. Do not make any changes to the medication of the person giving the sample without consulting their healthcare provider. If you are collecting a sample for a test that requires fasting, please consult with a healthcare provider, particularly if the person giving the sample has any medical conditions.
 - Do not pierce or collect blood from the smallest finger or any other body part not indicated in these instructions.
 - Do not use lancets on damaged, scarred, or infected skin.
 - Lancets may cause some bruising and tenderness at the puncture site. There is a small possibility of pain, nerve damage, scarring, bone damage, or infection. Consult a healthcare provider if you are concerned.
 - There is a risk of fainting with any blood sample. If you have a history of fainting or a medical condition that increases your risk of fainting, do not use this kit. Please consult your healthcare provider before attempting to collect your sample. If you proceed, ensure you are seated and have someone with you. If you feel faint or unwell during sample collection, stop and seek medical advice.
 - If assisting in child sample collection, the child should be seated and in a safe place while the lancet is being used. If the child has previously fainted while having blood collected or feels faint or unwell while their sample is being collected, please stop and consult with a healthcare provider.
 - Use only the contents included in this single-use kit for collecting the sample. Check the expiry date before you use the kit and do not use if expired.
 - If the kit or any of the contents are damaged or missing, the kit should not be used. Do not tamper with the kit contents — risk of serious injury.
 - Used kit contents can generally be safely disposed of in household waste. Refer to your local laws and guidelines.
 - Choking hazard - this kit contains small parts which may present a choking hazard to small children. Please ensure children are directly supervised while in the presence of the test kit and they do not put any contents near their face.
 - The sample must be returned in the kit box.
 - If you have difficulty using this kit due to vision impairment or difficulties in handling the test components, or understanding test instructions or results, please contact a healthcare provider.
- Do not use this kit if:**
- The person giving the sample has a bleeding disorder or is taking medication that increases the risk of bleeding.
 - The person giving the sample has impaired lymphatic drainage in both upper limbs.
 - The person giving the sample has a condition associated with poor healing, increased risk of infection, or skin ulceration — including but not limited to Raynaud's phenomenon, scleroderma, peripheral vascular disease, or peripheral neuropathy.
- Purpose/Intended use of device:**
The LetsGetChecked Blood Home Sample Collection Kit is an in vitro diagnostic medical device specifically intended for the primary containment and preservation of blood samples, for the purpose of in vitro diagnostic examination by our partner accredited laboratories.

PrivaPath Diagnostics Ltd.,
Unit 1, Northern Cross Business Park,
North Road, Dublin 11,
D11 XT26, Ireland.

Storage temperature: 15°C — 25°C.
Transport temperature: 4°C — 25°C.

For explanation of symbols see
[LetsGetChecked.com/symbol-glossary](https://letsgetchecked.com/symbol-glossary)

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PrivPath Diagnostics Ltd.,