LetsGetChecked Optum House Calls Urine  
Home Sample Collection IFU 240

SPECS AND INSTRUCTIONS

Colour Reference:

- LGC Teal: Pantone 7710 CP (CMYK)
- LGC Navy: Pantone 654 CP (CMYK)
- LGC Red: Pantone 1925 CP (CMYK)
- CMYK

Size:

108.5mm wide by 152.5mm high **folded**

Paper:

135gsm silk

Fold type:

4 page accordion

Font:

Proxima Nova  
(Regular, Semibold and Bold)

Part Number:

QR-IFU-240-R1 Dec-2024


\*Visual reference (showing the English side of the IFU) for how to fold down the sheet so the front page is the first one to see. Artwork and page count may vary.



FOLD GUIDE

----- Fold


Fold in a half first

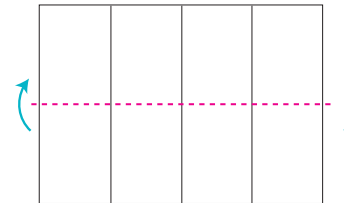
Then do the accordion fold with the text in English on the outside

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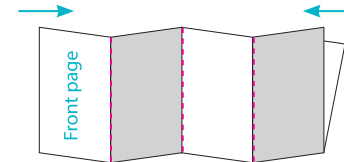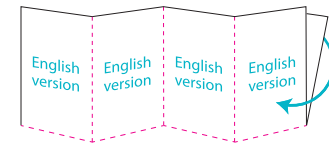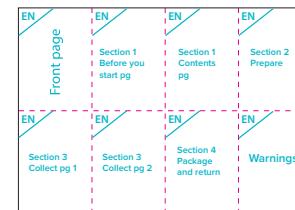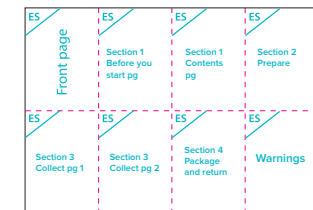QR-IFU-240-RT Dec 2024

## 1 Before you start

Consuming large amounts of protein before collecting your sample may impact your creatinine result by temporarily elevating it. This includes red meat, protein shakes, or supplements.

Do not collect a sample if you are on your period (menstruating).

Follow the instructions for filling the urine tube carefully. Tubes that are too full may be rejected.


Collect your sample in the morning and return on the same day.

Plan your sample return in advance. Check the return delivery card in your kit for details.

Scan the QR code below for an instructional video.
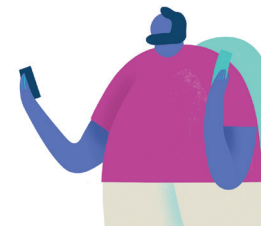
## Check your kit contents\*

Pipette

Urine sample tube

Urine collection cup

Biohazard bag


Return envelope

Kit box

\*The look or color of items may vary from what you see in these images.

## 2 Prepare for sample collection

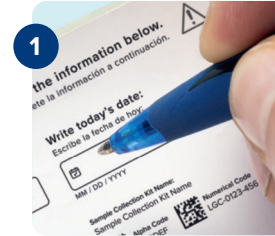1

Fill in the label on the biohazard bag.

2

Wash and dry your hands.

3

Remove all the contents from your kit box. Open the packaging.

↓ Flip over for next steps. ↓

## 3 Collect your urine sample

4

Open the collection cup by squeezing opposite corners. Do not tear or puncture the plastic lining.

5

Pass a small amount of your urine into the toilet and then pass your urine into the collection cup until it is half full.

6

Open the urine sample tube.

7

Use the pipette to slowly transfer your urine from the collection cup into the urine sample tube.

8

Slowly add urine to the tube. Stop when the tube is about half full. Remember, overfilled tubes may be rejected.

9

Screw the lid tightly onto the urine sample tube to close it. Check that the lid is tightly secured and does not leak. Wash and dry your hands.

10

Place the sample tube in the biohazard bag and seal the bag.

## 4 Package and return

11

Place the biohazard bag into the kit box. Double check you have filled in the biohazard bag label.

12

Place the kit box in the return envelope provided.

13

Send the kit back on the same day. Check the return delivery card in your kit box for your return method.

### WARNING:

- You must be over the age of 18 to use this home sample collection kit.
- The sample and information provided must be your own. Test results may not be released to a third party.
- Failure to follow these instructions may impact the accuracy of your test result or cause your sample to be rejected by the laboratory.
- Some medication and supplements may impact your test result. Do not make any changes to your medication without consulting your healthcare provider.
- Do not collect a sample if you are on your period (menstruating).
- Use only the contents included in this single-use kit for collecting your sample. Check the expiry date before you use the kit and do not use if expired.
- If the kit or any of the components are damaged or missing, the kit should not be used.
- Do not tamper with the kit contents – risk of serious injury.
- Used kit contents can generally be safely disposed of in household waste. Refer to your local laws and guidelines.

- Choking hazard – this kit contains small parts which may present a choking hazard. Keep away from children.
- The sample must be returned in the kit box.
- If you have difficulty using this kit due to vision impairment or difficulties in handling the test components, or understanding test instructions or results, please contact a healthcare provider.
- This test system is not a substitute for visits to a healthcare provider. The information provided by this product should not be used to start, stop, or change any course of treatment unless advised by your healthcare provider.

**Purpose/Intended use of device:**  
The LetsGetChecked Urine Home Sample Collection Kit is an in vitro diagnostic medical device specifically intended for the primary containment and preservation of urine samples, for the purpose of in vitro diagnostic examination by our partner accredited laboratories.

**Contact Customer Care:**  
1-855-554-1959 (TTY 711)

PrivaPath Diagnostics Ltd.,  
Unit 1, Northern Cross Business Park,  
North Road, Dublin 11,  
D11 XT26, Ireland.

Storage temperature: 15°C – 25°C.  
Transport temperature: 2°C – 30°C.

For explanation of symbols see  
[LetsGetChecked.com/symbol-glossary](https://www.letsgetchecked.com/symbol-glossary)

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